



Gasmonitor

Quality – Made in Germany

New Model



✓ CO₂

✓ N₂


Features:

- Gas monitoring for CO₂, N₂ and mixed gas (max. 21% O₂)
- Ensures continuous gas supply without interruption
- Controls and switches over automatically to standby cylinder if supply is exhausted
- Optical and acoustical alarm when switching over to standby cylinder



Features

- Measurement and monitoring of gas pressure
- Optical and acoustical alarm if pressure is too low or pressure is settled too high
- Automatic switch over to standby cylinder if supply is exhausted
- Graphic display with touch screen
- Quick couplers for gas connection
- Holder for installation at the wall

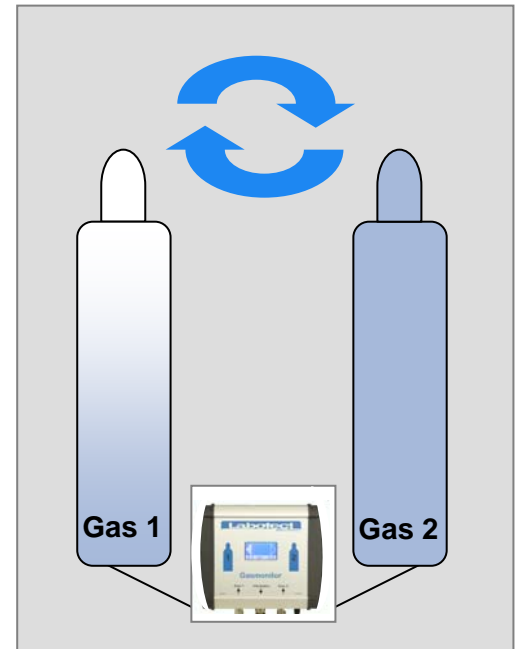

Details

- Pressure range: 0...2 bar
- Over pressure alarm: 1.5 bar
- Switch over at: 0.3 bar
- Remote alarm plug
- Power supply: 12 DC mains adapter plug

Scope of supply and Technical Data

- Gasmonitor
- Quick couplers and further connectors
- Holder for installation at the wall
- Mains adapter plug AC 100 - 240 V, DC 12 V 1.5 A
- Magnetic signs **GAS 1** and **GAS 2** for labelling of gas cylinder
- Dimensions: 19.5 x 5.3 x 18.1 cm (W x D x H)

Article no.: 14782



Specifications subject to change without notice

Labotect GmbH
Labor-Technik-Göttingen
Willi-Eichler-Str. 25
D-37079 Göttingen

Telefon +49 551 50501 0
Telefax +49 551 50501 11
e-mail info@labotect.com
Internet www.labotect.com

Amtsgericht Göttingen: HRB 1394
WEEE-Reg-Nr.: DE86039820
Geschäftsführerin / Managing Director
Delia Schinkel-Fleitmann

Rev 0_12/2009